

SURGICAL & GYNECOLOGY SUCTION UNIT The Best Just Got Better. DF-350A


- 1. The outer shall of the instrument adopted by ABS material. Clear simple lines and nice looking appearance create a new enjoyment for users. Humanized design with outstanding technology impressed users in big hospitals, medium/small clinics, nursing homes and home cares.
- 2. Adhering to the concept of environmental protection, ecological protection and econmical benefits makes each item of material be recycles as possible as it can to reduce waste and protect natural resource.

Functions:

- The Vacuum gauge is used to show the suction power, so users can easily read it.
- The suction tube directly connects to the suction bottle, so users can easily wash it. It's sanitary & safe.
- The handle attached to the suction bottle makes users easily lift and wash it.
- 360° Individual brake wheel features easy movement and high safety.
- The fastening rack installed at the side makes the suction tube be easily put.
- Surgical suction: Its smaller caliber is suitable for general vacuum.
- Obstetrical suction: Its bigger caliber is suitable for artificial abortion. Smaller caliber can be also used to suck baby out.
- Surgical & OBS. Select Valve on the panel can be switched, so there is no trouble of changing tube and it's more efficient.

The Instrument plane:

- Vacuum gauge
- Regulator valve
- O Power switch (On / Off)
- Surgical tube holder
- 6 Filter (L)
- Surgical suction adapter
- Surgical suction bottle (3000c.c.)
- Second (S)
 (S)
- Select valve (OBS. & Surgical)
- Obstetrical hose holder
- Floating control (L)
- Obstetrical suction adapter
- Obstetrical suction bottle (3000c.c.)

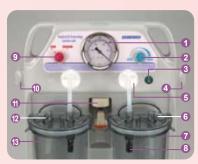





Filter (L)

Standard Accessories

1. Suction tube (S) 180cm 2 pc	3
2. Suction tube (L) 180cm ····· 1 pc	;
3. Suction catheter 1 pc	;
4. Foot switch 1 pc	;
5 Filter (L)	35

Specifications

250 W 1/3 HP Max. HP 10~720 mm/hg Vacuum Range Max. Flow 80 liter/min 3000c.c. + 3000c.c. Capacity **Machine Size** L45xW35xH84 (cm) L54xW44xH95 (cm) Packing Size Machine Weight 19.5 KGS **Gross weight** 25 KGS

The Side:

- 1 Power cord
- Adapter of foot switch
- 3 Circuit breaker
- 4 All casters with locking brakes
- 6 Foot switch

OIL-LESS VACUUM PUMP (90-026) Twin cylinder piston pump design.



Humanizing design: Low noise, Wear resistance.

Compact structure: Low power consumption, long life.

Fully oil-less: Non-pollution, free maintenance

Pump specification			
Voltage	120V	240V	
Current	2.4A	1.2A	
Capacitor	26uf	13uf	
R.P.M.	1720		
Poles	4		
Output	250W 1/3HP		
Max. flow rate	80 liter/min		
Max. suction	700~720 mm/hg		
Weight	9.0 kgs		
Thermal protector, auto reset			



1 Fix cover 2 Depositing box 3 Fan

