### **Accessories**

#### Standard









Safety Goggle

Fiber Cutter

Foot Switch

Fiber

### Optional



Dental Handpiece



Bleaching Handpiece



ENT Handpiece



Lipolysis Handpiece



ocus Handpiece



PM Handpiece



Fiber Stripper



Optitorque



Geide Wire







PLDD Sit





### Wuhan Gigaa Optronics Technology Co.,Ltd.

II-4-502,International Enterprise Center, Guanggu Road,Wuhan 430074,PRC Ph: 86-27-67848871 67848872 Fax: 86-27-67848873 Email: info@gigaalaser.com www.gigaalaser.com

# VFLAS TA



15W/30W/60W

Surgical Diode Laser Systems




www.gigaalaser.com


### The Innovative Platform Products

#### **VELAS Series Surgical Diode Laser Systems**

offer you a multiple choice for your different clinic requirement. Surgical applications require a technology that offers equally effective ablation and coagulation options. The VELAS series will offer you the wavelength options of 810nm, 940nm, 980nm and 1470nm, so that you could choose a laser which suits your need best.

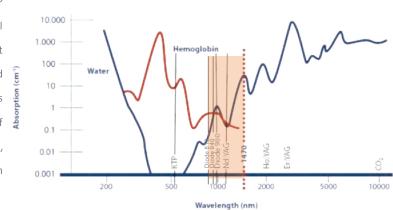

#### Lipolysis

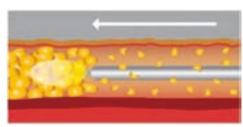

Laser lipolysis is the latest technology in fat reduction, which uses laser to dissolve fat, the cell membranes are broken up during this process. The thin fat lysate is released and it is either resorbed naturally or extracted manually by suction. The compatible operation kits for lipolysis by GIGAA LASER will provide a wide range of application in body and a totally new experience for patients.



#### PLDD(Percutaneous Laser Dsic Decompression)

Percutaneous laser disc decompression is the use of laser thermal effect and prominent vertebral Vaporization of disc nucleus and form part of the space, reducing the pressure inside the disc, from the Reduce or eliminates the oppression of the nerve, which achieve the purpose of treatment. The therapeutic method of the laser treatment for PLDD.



#### ENT(Ear, Nose and Throat)

- Dacryocystorhinostomy (DCR)
  - Turbinectomy
  - Nasal polypectomy
    - Tonsillectomy
  - Glossectomy
  - Endo sinus surgery
    - · Stapedectomy

### Technical Specifications

Laser Type	GaAIAs diode Laser		
Model	VELASII15A/B/D	VELASII30A/B	VELASII60A/B/C
Wavelength	810nm/980nm/1470nm	810nm/980nm	810nm/980nm/940nm
Maximum Power	15W	30W	60W
Operation Mode	CW, Single or Repeat Pulse		
Pulse Duration	10µs-3s		10ms-3s
Repetition Rate	0.2Hz-20KHz		0.2Hz-50Hz
Fibers	Fibers Of 200µm ,400µm And 600µm With SMA905 Connector		
	(200µm only could be used on VeLasII15B)		
Pilot Beam	Red Diode Laser of 635nm, Power<5mW		
Control Mode	True ColorTouch Screen		
Voltage/Current Rating	110/220 VAC,5A,50/60 Hz		
Dimensions	400(W)*385(L)*200(H) mm		
Weight	12.9Kg		13.4Kg
Safety Compliance	CE 0197		

## GIGAALASER ALWAYS PROVIDE A TOTAL SOLUTION TO YOU


The laser systems involve to Urology, EVLT, ENT, Dental, Gynecological, Orthopedic, Plastic surgery and other areas. Each year, we develop at least one new laser system, to bring the new technology and application, safety and health for the patients.

#### Dentistry

- · Internal decontamination
- · General surgery (cutting)
- ·Therapy
- · Tooth whitening
- · Surfaces of teeth
- · Internal gingivectomy

#### Phlebolog/Endovenous Laser Treatment

- · Great saphenous vein
- · Small saphenous vein
- · Perforating vein

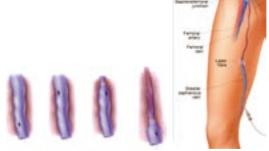



#### **Main Advantages**

Based on the experienced R.D, administration team, Gigaalaser Focuses on the medical diode laser system and accessories for application in a wide range of medical specialties. Using the perfect Chinese resource industry, we provide the high quality and the cost effective new type medical laser systems to the clients all over the world.



New friendly used protocol Big color touch screen

Output power calibration at the end of fibers



Highly integrated and modular design New design inside to make maintenance easy and low cost

High quality and reliability