

STERILE CONTAINER SYSTEMS PLASMA MODEL



1/1 FULL SIZE



3/4 QUARTER SIZE



1/2 HALF SIZE



DENTAL MODEL



Mini MODEL



1/2 DENTAL MODEL



ENDO MODEL



TRAY LAPAROSCOPY

Aygün plasma model container systems are designed for both steam and plasma sterilization (destroying the microorganisms by free radicals). The suitable accessories like paper filter for steam or polypropylene filter for plasma must be chosen by the user according to the need. The plasma model containers are equipped with carrying handles which have 90° mobility and its ergonomic design provides easy and secure carrying facilities. The gaskets inside the lid and filter retainer ensure perfect tightness while locking mechanism of this system provides excellent grip with its flexible structure.

Aygün plasma model sterilization container systems are manufactured with different length, width and height and can be supplied with perforated bottom and lids. All models and sizes have specific baskets and accessories which can be viewed in following parts.

STERILE CONTAINER SYSTEMS



LOW TEMPERATURE STERILIZATION CONTAINERS

For STERRAD[®]



* STERRAD is a registered trademark owned by Johnson & Johnson Company.

1/2 (HALF SIZE)



W: 280 mm

L: 310 mm

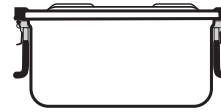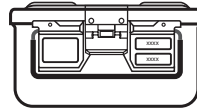


code	outside dimension	Specifications
K6-530Z	310 x 280 h:110 mm	Lid (Perforated)
K6-536Z	310 x 280 h:130 mm	
K6-542Z	310 x 280 h:150 mm	

Ps: All container measurements cover the outer size of the lids.

BOX (perforated)




outside dimension

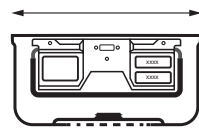

outside dimension



inside dimension

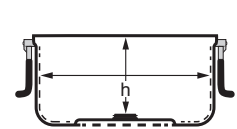

code	outside dimension	inside dimension	weight
K6-600Z	310 x 272 h:100 mm	253 x 258 h:83 mm	1385 g
K6-601Z	310 x 272 h:116 mm	253 x 258 h:99 mm	1417 g
K6-602Z	310 x 272 h:138 mm	253 x 258 h:121 mm	1590 g

PERFORATED LIDS



code	weight
K6-630Z	727 g


Dimension: 305 x 280 mm

SECURITY SEAL




K2-663Z
100 pcs. with indicator.
for hydrogen peroxide
plasma sterilization

POLYPROPYLENE FILTER

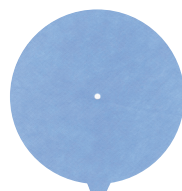


K2-741Z
Single use 100 pcs.
Diameter: 184 mm
for hydrogen peroxide
plasma sterilization

LABEL

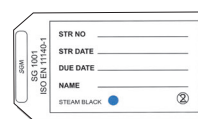

K2-672Z
100 pcs. with indicator.
for hydrogen peroxide
plasma sterilization

LABEL

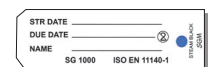

K2-673Z
100 pcs. with indicator.
for hydrogen peroxide
plasma sterilization

BASKET MADE FROM STAINLESS STEEL WIRE MESH

WIRE BASKET



code	outside dimension
K2-490	246 x 246 x h:30 mm
K2-491	246 x 246 x h:50 mm
K2-492	246 x 246 x h:70 mm
K2-493	246 x 246 x h:100 mm
K3-1045	ONLY LID, for baskets



WIRE BASKET WITH PERFORATED METAL PLATE

WIRE BASKET



code	outside dimension
K2-491 P	246 x 246 x h:50 mm
K2-492 P	246 x 246 x h:70 mm
K2-493 P	246 x 246 x h:100 mm
K3-1045 P	ONLY LID, for baskets




BASKET MADE FROM STAINLESS STEEL SHEET

BASKET (with stand)



	code	outside dimension
Box	K3-752	244 x 253 h:40 mm
	K3-753	244 x 253 h:60 mm
	K3-754	244 x 253 h:80 mm
	K3-755	244 x 253 h:110 mm
Lid	K3-763	ONLY LID, for baskets

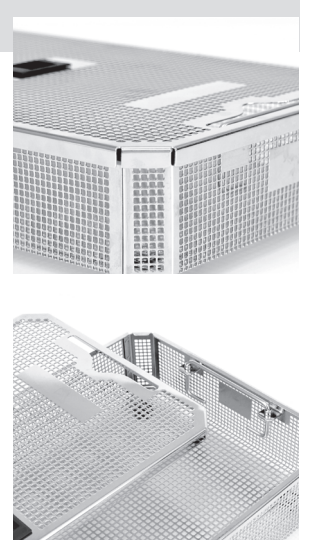


➤ The baskets which are without lid and less than 44 mm in height do not include the handles while the ones with lid include the handles.

STAND



code	height
K3-767	20 mm

STAND

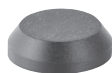

code	height
K3-766	5 mm



STANDARD
MODEL

29

BIO-BARRIER
MODEL

63

Polymer's PPSU
Lids

97

WINDOW
MODEL

119

B₃ Barrier
MODEL

139

PLASMA
MODEL

173

A1
MODEL

185

Dental, Mini,
Special

205

BASKETS

223

Autolock Mechanism
Emulating Indicator

243

ACCESSORIES

251