

STERILE CONTAINER SYSTEMS WINDOW MODEL



1/1 FULL SIZE



3/4 QUARTER SIZE



1/2 HALF SIZE

Aygün window model sterilization container system is one of its patented products. What makes this model different from the others is the window on the lid which its name comes from. This window is made of high performance PPSU transparent material which is resistant against high temperature, impacts, scratches, deformation, discoloration and chemicals and it opens users a way to visually check inside of the sterilization container without opening the lid. The window model lids assure great tightness due to its special assembly method and the gasket system in it.

As there is no filtering mechanism on the lid, window model sterilization container boxes strictly must be bottom perforated. Either single use paper or 2000 cycles PTFE Teflon Filters can be used in the box. In order to allow the air flow and reach to perfect steam penetration the filtering mechanism of this model must have enough space. For that purpose, this sterilization container model is supplied with a wire base inside the container as standard. Except for this container model, all Aygün boxes and lids are compatible with each other.

Aygün window model sterilization container systems are manufactured as five different length, width and height. All models and sizes have specific baskets and accessories which can be viewed in following parts. Modular auto lock mechanism (thermo lock) promising sterility assurance and indicator window providing visibility of the indicator card result by visual control without opening the lid are also available for this sterilization container model.

STERILE CONTAINER SYSTEMS



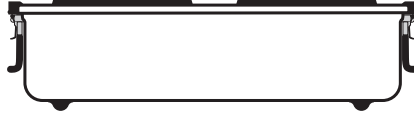
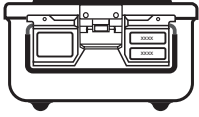
Window Model

» 1/1 Full Size


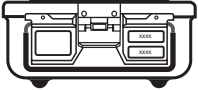







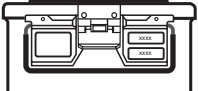






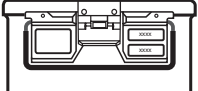













1/1 WINDOW MODEL

W: 285 mm

L: 600 mm




Shipped with wire base.

Specifications		Lid (Perforated)	
			
Height		Box (Perforated)	
		h:130 mm	
h:155 mm		K6-016D  K6-017D  K6-018D  K6-019D  K6-020D  K6-021D 	
h:175 mm		K6-022D  K6-023D  K6-024D  K6-025D  K6-026D  K6-027D 	
h:225 mm		K6-028D  K6-029D  K6-030D  K6-031D  K6-032D  K6-033D 	
h:275 mm		K6-034D  K6-035D  K6-036D  K6-037D  K6-038D  K6-039D 	



Ps: All container measurements cover the outer size of the lids.

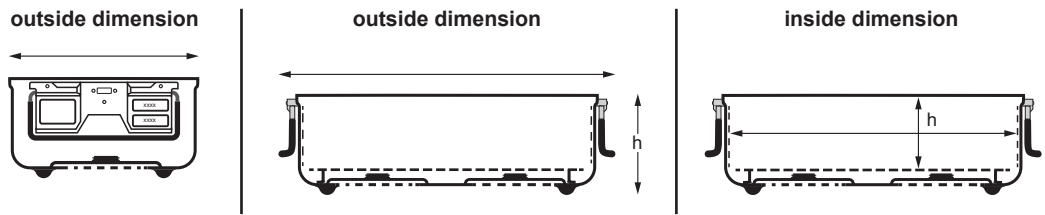BOX (perforated)



➡ Shipped with wire base.

code	outside dimension	inside dimension	weight
K6-100D	600 x 272 h:115 mm	543 x 258 h:66 mm	3003 g
K6-101D	600 x 272 h:131 mm	543 x 258 h:82 mm	3149 g
K6-102D	600 x 272 h:153 mm	543 x 258 h:104 mm	3307 g
K6-103D	600 x 272 h:205 mm	543 x 258 h:156 mm	3660 g
K6-104D	600 x 272 h:250 mm	543 x 258 h:201 mm	4077 g

LID WINDOW



Dimension: 602 x 285 mm

code	color	weight
K6-120D		1305 g
K6-121D		
K6-122D		
K6-123D		
K6-124D		
K6-125D		

STANDARD MODEL 29

BIO-BARRIER MODEL 63

Polymer's PPSU Lids 97

WINDOW MODEL 119

B₃ Barrier MODEL 139

PLASMA MODEL 173

A1 MODEL 185

Dental, Mini, Special 205


BASKETS 223

Autolock Mechanism Emulating Indicator 243


ACCESSORIES 251

BASKET MADE FROM STAINLESS STEEL WIRE MESH

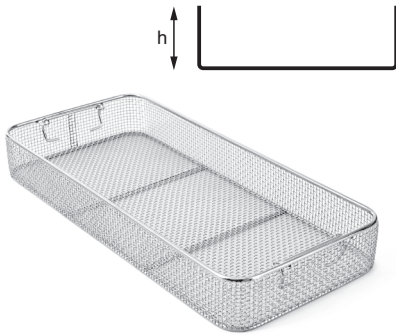
WIRE BASKET



code	outside dimension
K2-460	540 x 245 h:30 mm
K2-461	540 x 245 h:50 mm
K2-462	540 x 245 h:70 mm
K2-463	540 x 245 h:100 mm
K2-450	480 x 245 h:30 mm
K2-451	480 x 245 h:50 mm
K2-452	480 x 245 h:70 mm
K2-453	480 x 245 h:100 mm

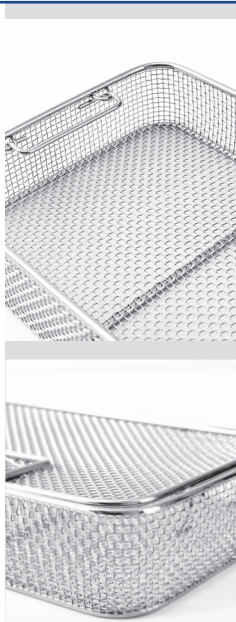

WIRE BASKET WITH LID

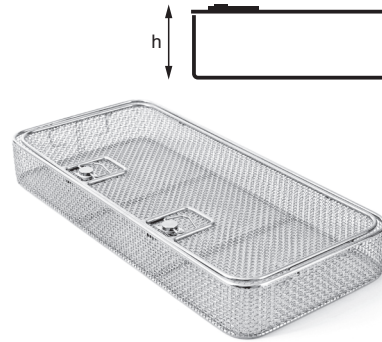

code	outside dimension
K3-1010	540 x 245 h:35 mm
K3-1011	540 x 245 h:55 mm
K3-1012	540 x 245 h:75 mm
K3-1013	540 x 245 h:105 mm
K3-1014	ONLY LID, for baskets
K3-1020	480 x 245 h:35 mm
K3-1021	480 x 245 h:55 mm
K3-1022	480 x 245 h:75 mm
K3-1023	480 x 245 h:105 mm
K3-1024	ONLY LID, for baskets




WIRE BASKET WITH PERFORATED METAL PLATE

WIRE BASKET



code	outside dimension
K2-461 P	540 x 245 h:50 mm
K2-462 P	540 x 245 h:70 mm
K2-463 P	540 x 245 h:100 mm
K2-451 P	480 x 245 h:50 mm
K2-452 P	480 x 245 h:70 mm
K2-453 P	480 x 245 h:100 mm

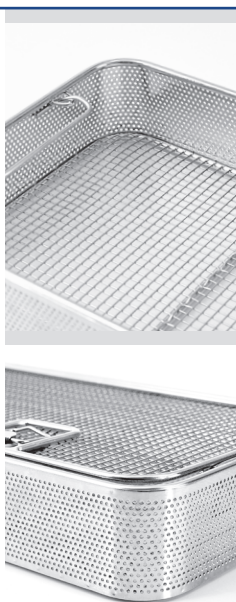

WIRE BASKET WITH LID







code	outside dimension
K3-1011 P	540 x 245 h:55 mm
K3-1012 P	540 x 245 h:75 mm
K3-1013 P	540 x 245 h:105 mm
K3-1014 P	ONLY LID, for baskets
K3-1021 P	480 x 245 h:55 mm
K3-1022 P	480 x 245 h:75 mm
K3-1023 P	480 x 245 h:105 mm
K3-1024 P	ONLY LID, for baskets





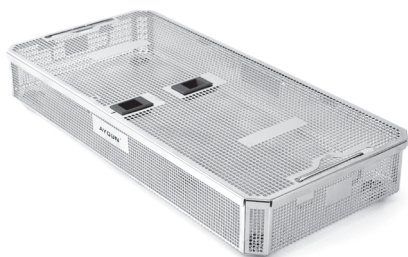BASKET MADE FROM STAINLESS STEEL SHEET


BASKET (with stand)



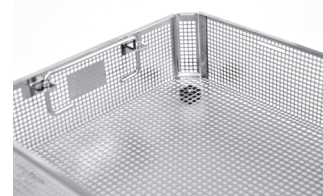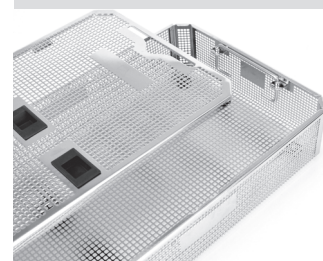


	code	outside dimension
Box and Lid 	K3-720	546 x 259 h:44 mm
	K3-721	546 x 259 h:64 mm
	K3-722	546 x 259 h:84 mm
	K3-723	546 x 259 h:114 mm
Box 	K3-740	540 x 253 h:40 mm
	K3-741	540 x 253 h:60 mm
	K3-742	540 x 253 h:80 mm
	K3-743	540 x 253 h:110 mm
	Lid 	K3-760

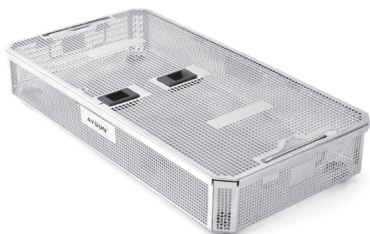


BASKET (with stand)



	code	outside dimension
Box and Lid 	K3-724	491 x 259 h:44 mm
	K3-725	491 x 259 h:64 mm
	K3-726	491 x 259 h:84 mm
	K3-727	491 x 259 h:114 mm
	Box 	K3-744
K3-745		485 x 253 h:60 mm
K3-746		485 x 253 h:80 mm
K3-747		485 x 253 h:110 mm
Lid 		K3-761


☞ The baskets which are without lid and less than 44 mm in height do not include the handles while the ones with lid include the handles.

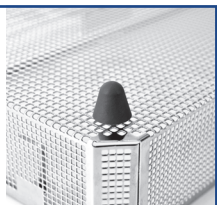STAND



code	height
K3-767	20 mm

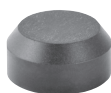STAND



code	height
K3-768	10 mm

STAND



code	height
K3-766	5 mm

IDENTIFICATION (aluminium)

LABEL HOLDER



code	color
K2-560	Red
K2-561	Blue
K2-562	Green
K2-563	Yellow
K2-564	Grey
K2-565	Black

LABELS



code	color
K2-570	Red
K2-571	Blue
K2-572	Green
K2-573	Yellow
K2-574	Grey
K2-575	Black

Coloured identification labels
Engraved to your specification



☞ Grey colored recommended for use in Sterrad.
Color fading By effect of Hydrogen Peroxide Gas Plasma on other colors.

STANDARD MODEL	29
BIO-BARRIER MODEL	63
Polymer's PPSU Lids	97
WINDOW MODEL	119
B ₃ Barrier MODEL	139
PLASMA MODEL	173
A1 MODEL	185
Dental, Mini, Special	205
BASKETS	223
Autolock Mechanism Emulating Indicator	243
ACCESSORIES	251

DIVIDERS



K2-428
540 x 250 x 130 mm
K2-429
540 x 250 x 180 mm

WIRE BASE



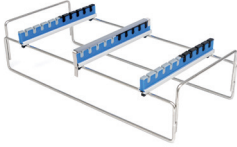
K2-440
545 x 255 x 30 mm

DIVIDERS




K2-425
540 x 250 x 130 mm
K2-426
540 x 250 x 180 mm

RACK (for endoscopy)



K3-1090
505 x 250 x 125 mm
• Holds 10 endoscopic instruments
• Fits into a 1/1 size 20 cm deep container
• Wire basket 1/1 50 mm
• Container 1/1 200 mm
• Silicone mat 520 x 230 mm

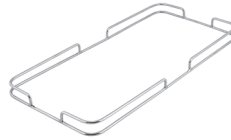DRAPE RETAINER



K2-420
545 x 255 x 30 mm

DRAPE RETAINER



K2-400
545 x 255 mm

ALUMINIUM BASKET (with lid)



K3-1119
245 x 130 x 55 mm

ALUMINIUM BASKET (with lid)



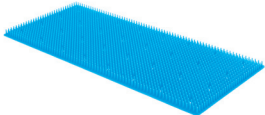
K3-1120
245 x 130 x 70 mm

ALUMINIUM BASKET (with lid)




K3-1121
245 x 130 x 105 mm

SILICONE MAT



K2-526
520 x 230 mm
K2-528
475 x 230 mm

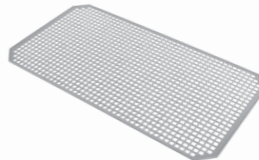SILICONE MESH MAT



K2-531
485 x 250 mm
K2-533
534 x 247 mm

SILICONE HANDLE



K2-710
K2-711
K2-712
K2-713
K2-714
K2-715

PAPER FILTER



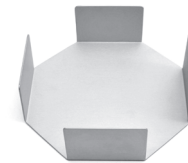
K2-740
Single use 100 pcs with indicator. Aygun mark.
• Diameter: 184 mm

TEFLON FILTER



K2-742
Number of sterilization cycle > 2000
• Diameter: 194 mm

FILTER DISPENSER




K2-730

IDENTIFICATION LABELS (aluminium)




K2-700
K2-701
K2-702
K2-703
K2-704
K2-705

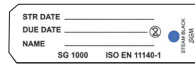LABELS (paper)



K2-673
100 pcs. With indicator.
For small size and special containers

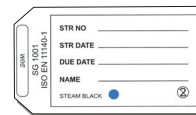LABELS (paper)



K2-672
100 pcs. With indicator.
For full and Half size containers

(-)

AUTO LOCK MECHANISM



K1-900
(2 Pcs)

EMULATING INDICATOR CARD (class 6)



K2-674
Paper 250 Pcs.



SECURITY SEAL (blue)
K2-663
100 pcs. with indicator. made of plastic



SECURITY SEAL (blue)
K2-664
100 pcs. made of plastic



SECURITY SEAL (blue)
K2-661
100 pcs. with number. made of plastic



SECURITY SEAL (red)
K2-662
100 pcs. with indicator. made of plastic



SECURITY SEAL (red)
K2-660
100 pcs. made of plastic



SECURITY SEAL (red)
K2-665
100 pcs. with number. made of plastic



DISPOSAL STAMP
K2-717
40 x 35 x 9 mm
Requested number or letter must be informed during the order.



LABEL HOLDER
K2-716
50 x 35 x 13 mm
Holder for identification labels

(-)