

STERILE CONTAINER SYSTEMS STANDARD MODEL



1/1 FULL SIZE



3/4 QUARTER SIZE



1/2 HALF SIZE



EXTRA LARGE SIZE



EXTRA LONG SIZE

Aygün standard model container systems are designed for perfect steam flow and penetration for loads suitable to the standards and can be used with single use paper filter or reusable PTFE Teflon filters. The standard model sterilization containers are equipped with silicone coated carrying handles which have 90° mobility and its ergonomic design provides easy and secure carrying facilities. The gaskets inside the lid and filter retainer ensure perfect tightness while locking mechanism of this system provides excellent grip with its flexible structure.

Aygün standard model sterilization container systems are manufactured as five different length, width and height and can be supplied with perforated or non-perforated boxes, perforated lids and safety lids for effective external protection. These safety lids do not interfere with sterilization and they provides protection to the perforated lid and the filter system against hard objects, liquids and dust as a second lid. All models and sizes have specific baskets and accessories which can be viewed in following parts. Modular auto lock mechanism (thermo lock) promising sterility assurance and indicator window providing visibility of the indicator card result by visual control without opening the lid are also available for this sterilization container model.

STERILE CONTAINER SYSTEMS





Standard Model

» *1/2 Half Size*

1/2 STANDARD MODEL

W: 280 mm

L: 310 mm

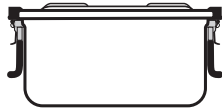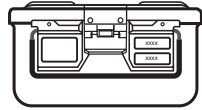


Specifications		Lid (Perforated)		Safety Lid	
		Box (Non Perforated)	Box (Perforated)	Box (Non Perforated)	Box (Perforated)
H: 110 mm		K6-500	K6-530	K6-640	K6-670
		K6-501	K6-531	K6-641	K6-671
		K6-502	K6-532	K6-642	K6-672
		K6-503	K6-533	K6-643	K6-673
		K6-504	K6-534	K6-644	K6-674
		K6-505	K6-535	K6-645	K6-675
H: 130 mm		K6-506	K6-536	K6-646	K6-676
		K6-507	K6-537	K6-647	K6-677
		K6-508	K6-538	K6-648	K6-678
		K6-509	K6-539	K6-649	K6-679
		K6-510	K6-540	K6-650	K6-680
		K6-511	K6-541	K6-651	K6-681
H: 150 mm		K6-512	K6-542	K6-652	K6-682
		K6-513	K6-543	K6-653	K6-683
		K6-514	K6-544	K6-654	K6-684
		K6-515	K6-545	K6-655	K6-685
		K6-516	K6-546	K6-656	K6-686
		K6-517	K6-547	K6-657	K6-687
H: 200 mm		K6-518	K6-548	K6-658	K6-688
		K6-519	K6-549	K6-659	K6-689
		K6-520	K6-550	K6-660	K6-690
		K6-521	K6-551	K6-661	K6-691
		K6-522	K6-552	K6-662	K6-692
		K6-523	K6-553	K6-663	K6-693
H: 250 mm		K6-524	K6-554	K6-664	K6-694
		K6-525	K6-555	K6-665	K6-695
		K6-526	K6-556	K6-666	K6-696
		K6-527	K6-557	K6-667	K6-697
		K6-528	K6-558	K6-668	K6-698
		K6-529	K6-559	K6-669	K6-699

Ps: All container measurements cover the outer size of the lids.


BOX (perforated)




outside dimension

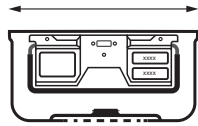

outside dimension

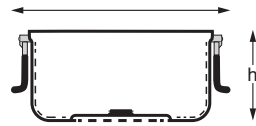

inside dimension

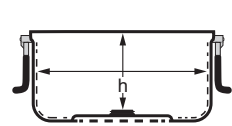


code	outside dimension	inside dimension	weight
K6-600	310 x 272 h:100 mm	253 x 258 h:83 mm	1385 g
K6-601	310 x 272 h:116 mm	253 x 258 h:99 mm	1417 g
K6-602	310 x 272 h:138 mm	253 x 258 h:121 mm	1590 g
K6-603	310 x 272 h:190 mm	253 x 258 h:173 mm	1821 g
K6-604	310 x 272 h:235 mm	253 x 258 h:218 mm	1947 g


BOX (non perforated)




outside dimension

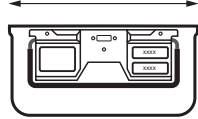

outside dimension

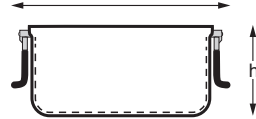

inside dimension

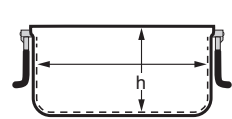

code	outside dimension	inside dimension	weight
K6-605	310 x 272 h:100 mm	253 x 258 h:94 mm	1267 g
K6-606	310 x 272 h:116 mm	253 x 258 h:110 mm	1315 g
K6-607	310 x 272 h:138 mm	253 x 258 h:132 mm	1427 g
K6-608	310 x 272 h:190 mm	253 x 258 h:184 mm	1660 g
K6-609	310 x 272 h:235 mm	253 x 258 h:229 mm	1940 g

PERFORATED LIDS



code	color	weight
K6-630		727 g
K6-631		
K6-632		
K6-633		
K6-634		
K6-635		

Dimension: 312 x 280 mm

REMOVAL LIDS (non perforated)



code	color	weight
K6-620		629 g
K6-621		
K6-622		
K6-623		
K6-624		
K6-625		

Dimension: 312 x 280 mm

SAFETY LIDS



code	color	weight
K6-700		775 g
K6-701		
K6-702		
K6-703		
K6-704		
K6-705		

Dimension: 317 x 285 mm

STANDARD MODEL

29

BIO-BARRIER MODEL

63

Polymer's PPSU Lids

97

WINDOW MODEL

119

B₃ Barrier MODEL

139

PLASMA MODEL

173

A1 MODEL

185

Dental, Mini, Special

205

BASKETS

223

Autolock Mechanism Emulating Indicator


243

ACCESSORIES

251


BASKET MADE FROM STAINLESS STEEL WIRE MESH

WIRE BASKET

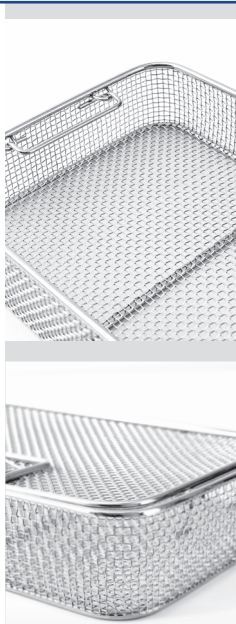

code	outside dimension
K2-490	246 x 246 x h:30 mm
K2-491	246 x 246 x h:50 mm
K2-492	246 x 246 x h:70 mm
K2-493	246 x 246 x h:100 mm


WIRE BASKET WITH LID

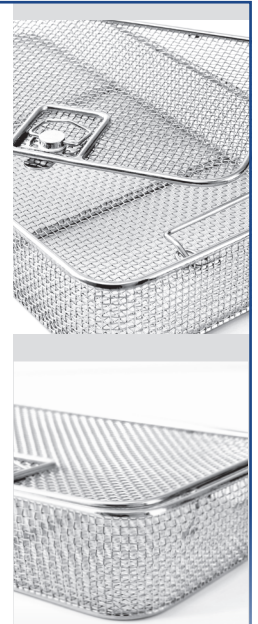

code	outside dimension
K3-1040	246 x 246 h:35 mm
K3-1041	246 x 246 h:55 mm
K3-1042	246 x 246 h:75 mm
K3-1043	246 x 246 h:105 mm
K3-1045	ONLY LID, for baskets

WIRE BASKET WITH PERFORATED METAL PLATE

WIRE BASKET

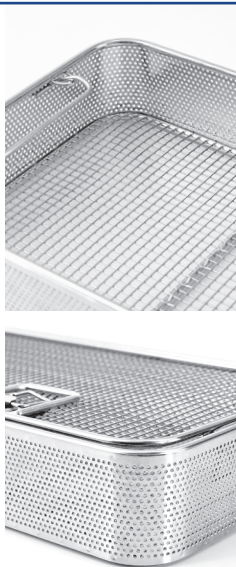

code	outside dimension
K2-491 P	246 x 246 x h:50 mm
K2-492 P	246 x 246 x h:70 mm
K2-493 P	246 x 246 x h:100 mm

WIRE BASKET WITH LID









code	outside dimension
K3-1041 P	246 x 246 x h:55 mm
K3-1042 P	246 x 246 x h:75 mm
K3-1043 P	246 x 246 x h:105 mm
K3-1045 P	ONLY LID, for baskets

BASKET MADE FROM STAINLESS STEEL SHEET

BASKET (with stand)

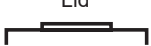

	code	outside dimension
Box and Lid 	K3-732	250 x 259 h:44 mm
	K3-733	250 x 259 h:64 mm
	K3-734	250 x 259 h:84 mm
	K3-735	250 x 259 h:114 mm
Box 	K3-752	244 x 253 h:40 mm
	K3-753	244 x 253 h:60 mm
	K3-754	244 x 253 h:80 mm
	K3-755	244 x 253 h:110 mm
Lid 	K3-763	ONLY LID, for baskets

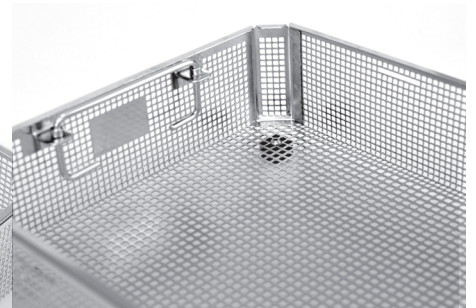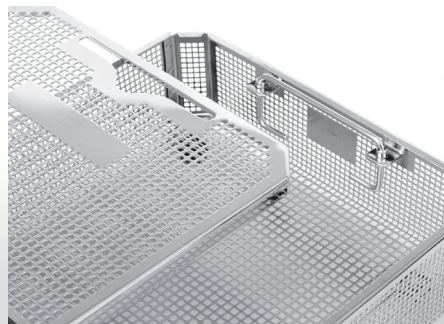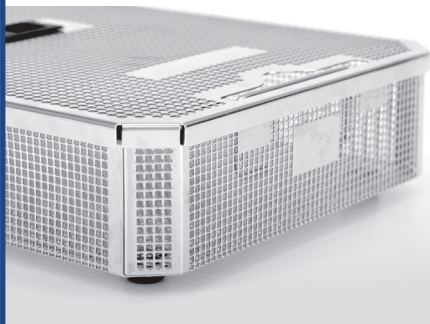


↻ The baskets which are without lid and less than 44 mm in height do not include the handles while the ones with lid include the handles.

STAND





code	height
K3-767	20 mm

STAND

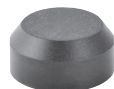

code	height
K3-768	10 mm

STAND

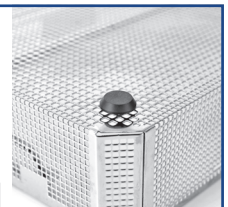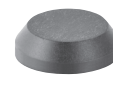

code	height
K3-766	5 mm

IDENTIFICATION (aluminium)

LABEL HOLDER



code	color
K2-560	Red
K2-561	Blue
K2-562	Green
K2-563	Yellow
K2-564	Grey
K2-565	Black

LABELS



code	color
K2-570	Red
K2-571	Blue
K2-572	Green
K2-573	Yellow
K2-574	Grey
K2-575	Black


Coloured identification labels
Engraved to your specification



↻ Grey colored recommended for use in Sterrad.
Color fading By effect of Hydrogen Peroxide Gas Plasma on other colors.

STANDARD MODEL	29
BIO-BARRIER MODEL	63
Polymer's PPSU Lids	97
WINDOW MODEL	119
B ₃ Barrier MODEL	139
PLASMA MODEL	173
A1 MODEL	185
Dental, Mini, Special	205
BASKETS	223
Autolock Mechanism Emulating Indicator	243
ACCESSORIES	251

DIVIDERS

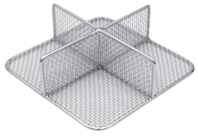

K2-483
255 x 245 x 130 mm
K2-484
255 x 245 x 180 mm

WIRE BASE



K2-482
255 x 255 x 30 mm


(-)

RACK (for kerrison)



K3-1094
220 x 250 x 160 mm
• Holds 10 kerrison punches
• Fits into a 1/2 and 1/1 size
• 20 cm deep container
• Wire basket 1/2 50 mm
• Wire basket 1/1 50 mm
• Container 1/2 200 mm
• Container 1/1 200 mm

DRAPE RETAINER

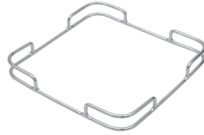

K3-481
255 x 255 x 30 mm

DRAPE RETAINER



K2-480
255 x 255 mm

ALUMINIUM BASKET (with lid)



K3-1119
245 x 130 x 55 mm

ALUMINIUM BASKET (with lid)



K3-1120
245 x 130 x 70 mm

ALUMINIUM BASKET (with lid)




K3-1121
245 x 130 x 105 mm

SILICONE MAT



K2-521
250 x 235 mm

SILICONE MESH MAT

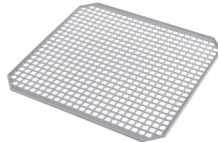

K2-530
250 x 240 mm

SILICONE HANDLE



K2-710
K2-711
K2-712
K2-713
K2-714
K2-715

PAPER FILTER




K2-740
Single use 100 pcs with indicator. Aygun mark.
• Diameter: 184 mm

TEFLON FILTER



K2-742
Number of sterilization cycle > 2000
• Diameter: 194 mm

FILTER DISPENSER

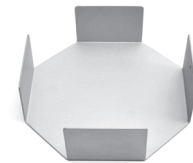


K2-730

IDENTIFICATION LABELS (aluminium)




K2-700
K2-701
K2-702
K2-703
K2-704
K2-705

LABELS (paper)

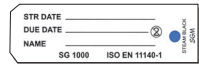

K2-673
100 pcs. With indicator.
For small size and special containers

LABELS (paper)

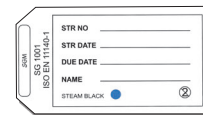

K2-672
100 pcs. With indicator.
For full and Half size containers


(-)

AUTO LOCK MECHANISM



K1-900
(2 Pcs)

EMULATING INDICATOR CARD (class 6)

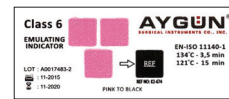

K2-674
Paper 250 Pcs.



SECURITY SEAL (blue)
K2-663
100 pcs. with indicator. made of plastic



SECURITY SEAL (blue)
K2-664
100 pcs. made of plastic



SECURITY SEAL (blue)
K2-661
100 pcs. with number. made of plastic



SECURITY SEAL (red)
K2-662
100 pcs. with indicator. made of plastic



SECURITY SEAL (red)
K2-660
100 pcs. made of plastic



SECURITY SEAL (red)
K2-665
100 pcs. with number. made of plastic



DISPOSAL STAMP
K2-717
40 x 35 x 9 mm
Requested number or letter must be informed during the order.



LABEL HOLDER
K2-716
50 x 35 x 13 mm
Holder for identification labels

(-)

Assembly and disassembly method of safety lid for safe model containers.

