

STERILE CONTAINER SYSTEMS WINDOW MODEL



1/1 FULL SIZE



3/4 QUARTER SIZE



1/2 HALF SIZE

Aygün window model sterilization container system is one of its patented products. What makes this model different from the others is the window on the lid which its name comes from. This window is made of high performance PPSU transparent material which is resistant against high temperature, impacts, scratches, deformation, discoloration and chemicals and it opens users a way to visually check inside of the sterilization container without opening the lid. The window model lids assure great tightness due to its special assembly method and the gasket system in it.

As there is no filtering mechanism on the lid, window model sterilization container boxes strictly must be bottom perforated. Either single use paper or 2000 cycles PTFE Teflon Filters can be used in the box. In order to allow the air flow and reach to perfect steam penetration the filtering mechanism of this model must have enough space. For that purpose, this sterilization container model is supplied with a wire base inside the container as standard. Except for this container model, all Aygün boxes and lids are compatible with each other.

Aygün window model sterilization container systems are manufactured as five different length, width and height. All models and sizes have specific baskets and accessories which can be viewed in following parts. Modular auto lock mechanism (thermo lock) promising sterility assurance and indicator window providing visibility of the indicator card result by visual control without opening the lid are also available for this sterilization container model.

STERILE CONTAINER SYSTEMS



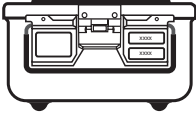
Window Model

» *1/2 Half Size*


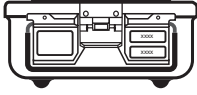
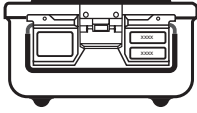
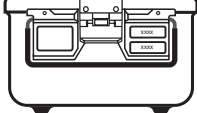
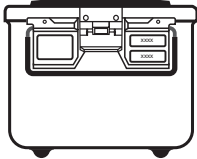
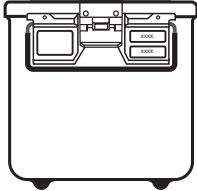
1/2 WINDOW MODEL

W: 280 mm

L: 310 mm



➡ Shipped with wire base.


Specifications ➡		Lid (Perforated)
Height		
		Box (Perforated)
h:130 mm		K6-500D <input type="checkbox"/>
		K6-501D <input type="checkbox"/>
		K6-502D <input type="checkbox"/>
		K6-503D <input type="checkbox"/>
		K6-504D <input type="checkbox"/>
		K6-505D <input type="checkbox"/>
h:150 mm		K6-506D <input type="checkbox"/>
		K6-507D <input type="checkbox"/>
		K6-508D <input type="checkbox"/>
		K6-509D <input type="checkbox"/>
		K6-510D <input type="checkbox"/>
h:170 mm		K6-511D <input type="checkbox"/>
		K6-512D <input type="checkbox"/>
		K6-513D <input type="checkbox"/>
		K6-514D <input type="checkbox"/>
		K6-515D <input type="checkbox"/>
h:220 mm		K6-516D <input type="checkbox"/>
		K6-517D <input type="checkbox"/>
		K6-518D <input type="checkbox"/>
		K6-519D <input type="checkbox"/>
		K6-520D <input type="checkbox"/>
h:270 mm		K6-521D <input type="checkbox"/>
		K6-522D <input type="checkbox"/>
		K6-523D <input type="checkbox"/>
		K6-524D <input type="checkbox"/>
		K6-525D <input type="checkbox"/>
		K6-526D <input type="checkbox"/>
		K6-527D <input type="checkbox"/>
		K6-528D <input type="checkbox"/>
		K6-529D <input type="checkbox"/>




➡ Ps: All container measurements cover the outer size of the lids.

BOX (perforated)

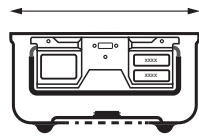
outside dimension



outside dimension

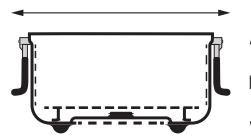

inside dimension

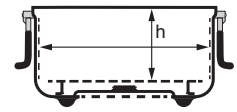

➔ Shipped with wire base.

code	outside dimension	inside dimension	weight
K6-600D	310 x 272 h:115 mm	253 x 258 h:66 mm	1801 g
K6-601D	310 x 272 h:131 mm	253 x 258 h:82 mm	1833 g
K6-602D	310 x 272 h:153 mm	253 x 258 h:104 mm	2006 g
K6-603D	310 x 272 h:205 mm	253 x 258 h:156 mm	2237 g
K6-604D	310 x 272 h:250 mm	253 x 258 h:201 mm	2363 g

LID WINDOW



Dimension: 312 x 280 mm

code	color	weight
K6-630D		736 g
K6-631D		
K6-632D		
K6-633D		
K6-634D		
K6-635D		

STANDARD MODEL 29

BIO-BARRIER MODEL 63

Polymer's PPSU Lids 97

WINDOW MODEL 119

B₃ Barrier MODEL 139

PLASMA MODEL 173

A1 MODEL 185

Dental, Mini, Special 205


BASKETS 223

Autolock Mechanism Emulating Indicator 243

ACCESSORIES 251

BASKET MADE FROM STAINLESS STEEL WIRE MESH

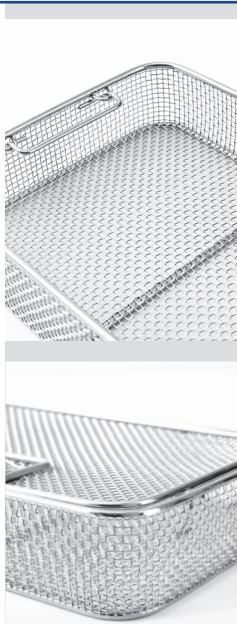WIRE BASKET



code	outside dimension
K2-490	246 x 246 x h:30 mm
K2-491	246 x 246 x h:50 mm
K2-492	246 x 246 x h:70 mm
K2-493	246 x 246 x h:100 mm


WIRE BASKET WITH LID



code	outside dimension
K3-1040	246 x 246 h:35 mm
K3-1041	246 x 246 h:55 mm
K3-1042	246 x 246 h:75 mm
K3-1043	246 x 246 h:105 mm
K3-1045	ONLY LID, for baskets

WIRE BASKET WITH PERFORATED METAL PLATE

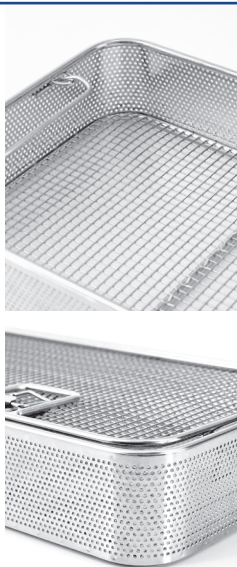WIRE BASKET



code	outside dimension
K2-491 P	246 x 246 x h:50 mm
K2-492 P	246 x 246 x h:70 mm
K2-493 P	246 x 246 x h:100 mm

WIRE BASKET WITH LID









code	outside dimension
K3-1041 P	246 x 246 x h:55 mm
K3-1042 P	246 x 246 x h:75 mm
K3-1043 P	246 x 246 x h:105 mm
K3-1045 P	ONLY LID, for baskets

BASKET MADE FROM STAINLESS STEEL SHEET

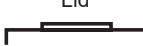BASKET (with stand)




	code	outside dimension
Box and Lid 	K3-732	250 x 259 h:44 mm
	K3-733	250 x 259 h:64 mm
	K3-734	250 x 259 h:84 mm
	K3-735	250 x 259 h:114 mm
Box 	K3-752	244 x 253 h:40 mm
	K3-753	244 x 253 h:60 mm
	K3-754	244 x 253 h:80 mm
	K3-755	244 x 253 h:110 mm
Lid 	K3-763	ONLY LID, for baskets

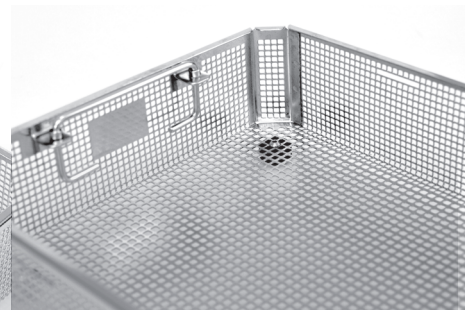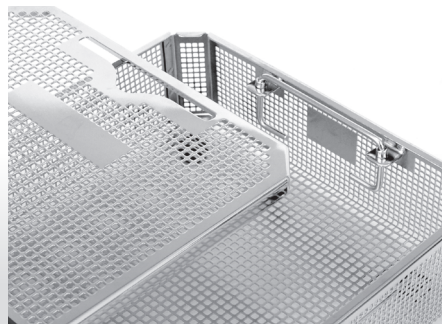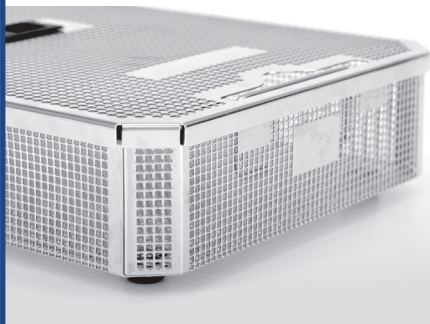

➤ The baskets which are without lid and less than 44 mm in height do not include the handles while the ones with lid include the handles.

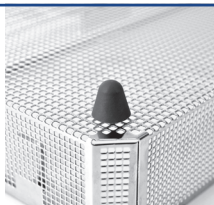STAND



code	height
K3-767	20 mm

STAND



code	height
K3-768	10 mm

STAND



code	height
K3-766	5 mm

IDENTIFICATION (aluminium)

LABEL HOLDER



code	color
K2-560	Red
K2-561	Blue
K2-562	Green
K2-563	Yellow
K2-564	Grey
K2-565	Black

LABELS



code	color
K2-570	Red
K2-571	Blue
K2-572	Green
K2-573	Yellow
K2-574	Grey
K2-575	Black

Coloured identification labels
Engraved to your specification



➤ Grey colored recommended for use in Sterrad.
Color fading By effect of Hydrogen Peroxide Gas Plasma on other colors.

STANDARD MODEL

29

BIO-BARRIER MODEL

63

Polymer's PPSU Lids

97

WINDOW MODEL

119

B₃ Barrier MODEL

139

PLASMA MODEL

173

A1 MODEL

185

Dental, Mini, Special

205

BASKETS

223

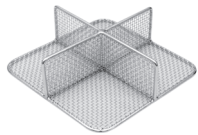
Autolock Mechanism Emulating Indicator

243

ACCESSORIES

251

DIVIDERS



K2-483
255 x 245 x 130 mm
K2-484
255 x 245 x 180 mm

WIRE BASE



K2-482
255 x 255 x 30 mm

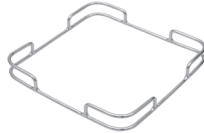
(-)

RACK (for kerrison)



K3-1094
220 x 250 x 160 mm
• Holds 10 kerrison punches
• Fits into a 1/2 and 1/1 size
• 20 cm deep container
• Wire basket 1/2 50 mm
• Wire basket 1/1 50 mm
• Container 1/2 200 mm
• Container 1/1 200 mm

DRAPE RETAINER



K3-481
255 x 255 x 30 mm

DRAPE RETAINER



K2-480
255 x 255 mm

ALUMINIUM BASKET (with lid)



K3-1119
245 x 130 x 55 mm

ALUMINIUM BASKET (with lid)



K3-1120
245 x 130 x 70 mm

ALUMINIUM BASKET (with lid)



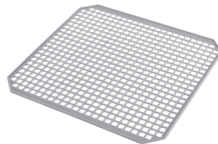
K3-1121
245 x 130 x 105 mm

SILICONE MAT



K2-521
250 x 235 mm







SILICONE MESH MAT



K2-530
250 x 240 mm

SILICONE HANDLE



K2-710 
K2-711 
K2-712 
K2-713 
K2-714 
K2-715 

PAPER FILTER



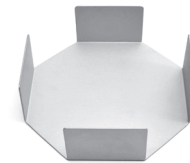
K2-740
Single use 100 pcs with indicator. Aygun mark.
• Diameter: 184 mm

TEFLON FILTER



K2-742
Number of sterilization cycle > 2000
• Diameter: 194 mm







FILTER DISPENSER



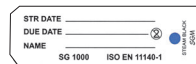
K2-730

IDENTIFICATION LABELS (aluminium)




K2-700 
K2-701 
K2-702 
K2-703 
K2-704 
K2-705 

LABELS (paper)



K2-673
100 pcs. With indicator.
For small size and special containers

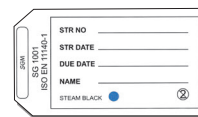LABELS (paper)



K2-672
100 pcs. With indicator.
For full and Half size containers

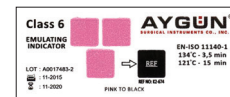
(-)

AUTO LOCK MECHANISM



K1-900
(2 Pcs)

EMULATING INDICATOR CARD (class 6)



K2-674
Paper 250 Pcs.



SECURITY SEAL (blue)
K2-663
100 pcs. with indicator. made of plastic



SECURITY SEAL (blue)
K2-664
100 pcs. made of plastic



SECURITY SEAL (blue)
K2-661
100 pcs. with number. made of plastic



SECURITY SEAL (red)
K2-662
100 pcs. with indicator. made of plastic



SECURITY SEAL (red)
K2-660
100 pcs. made of plastic



SECURITY SEAL (red)
K2-665
100 pcs. with number. made of plastic



DISPOSAL STAMP
K2-717
40 x 35 x 9 mm
Requested number or letter must be informed during the order.



LABEL HOLDER
K2-716
50 x 35 x 13 mm
Holder for identification labels

(-)