

## STERILE CONTAINER SYSTEMS WINDOW MODEL



1/1 FULL SIZE



3/4 QUARTER SIZE



1/2 HALF SIZE

Aygün window model sterilization container system is one of its patented products. What makes this model different from the others is the window on the lid which its name comes from. This window is made of high performance PPSU transparent material which is resistant against high temperature, impacts, scratches, deformation, discoloration and chemicals and it opens users a way to visually check inside of the sterilization container without opening the lid. The window model lids assure great tightness due to its special assembly method and the gasket system in it.

As there is no filtering mechanism on the lid, window model sterilization container boxes strictly must be bottom perforated. Either single use paper or 2000 cycles PTFE Teflon Filters can be used in the box. In order to allow the air flow and reach to perfect steam penetration the filtering mechanism of this model must have enough space. For that purpose, this sterilization container model is supplied with a wire base inside the container as standard. Except for this container model, all Aygün boxes and lids are compatible with each other.

Aygün window model sterilization container systems are manufactured as five different length, width and height. All models and sizes have specific baskets and accessories which can be viewed in following parts. Modular auto lock mechanism (thermo lock) promising sterility assurance and indicator window providing visibility of the indicator card result by visual control without opening the lid are also available for this sterilization container model.

**STERILE CONTAINER SYSTEMS**





**Window Model**

**» 3/4 Quarter Size**

### 3/4 WINDOW MODEL

W: 285 mm


L: 475 mm

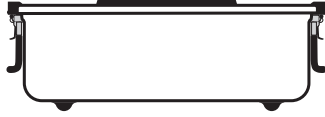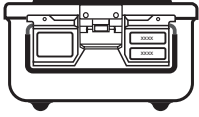

Specifications		Lid (Perforated)
		Box (Perforated)
h: 130 mm		<b>K6-240D</b>
		<b>K6-241D</b>
		<b>K6-242D</b>
		<b>K6-243D</b>
		<b>K6-244D</b>
h: 150 mm		<b>K6-245D</b>
		<b>K6-246D</b>
		<b>K6-247D</b>
		<b>K6-248D</b>
		<b>K6-249D</b>
h: 175 mm		<b>K6-250D</b>
		<b>K6-251D</b>
		<b>K6-252D</b>
		<b>K6-253D</b>
		<b>K6-254D</b>
		<b>K6-255D</b>
		<b>K6-256D</b>
		<b>K6-257D</b>




Shipped with wire base.

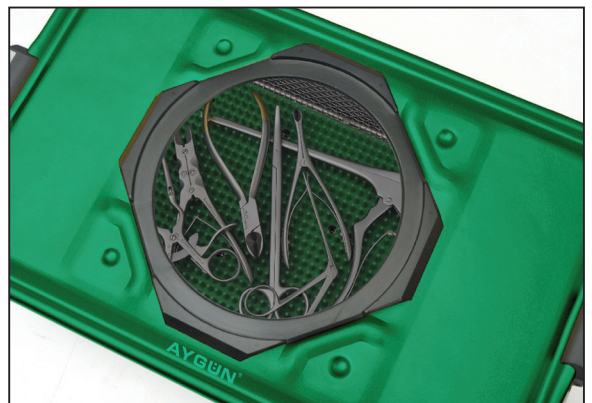


Ps: All container measurements cover the outer size of the lids.

**BOX** (perforated)


outside dimension

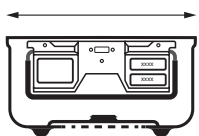

outside dimension

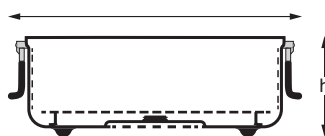

inside dimension

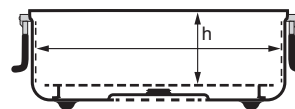

Shipped with wire base.

code	outside dimension	inside dimension	weight
<b>K6-300D</b>	475 x 272 h:115 mm	422 x 258 h:66 mm	2363 g
<b>K6-301D</b>	475 x 272 h:131 mm	422 x 258 h:82 mm	2407 g
<b>K6-302D</b>	475 x 272 h:153 mm	422 x 258 h:104 mm	2551 g
<b>K6-303D</b>	475 x 272 h:200 mm	422 x 258 h:151 mm	2811 g

**LID WINDOW**




Dimension: 477 x 285 mm

code	color	weight
<b>K6-340D</b>		918 g
<b>K6-341D</b>		
<b>K6-342D</b>		
<b>K6-343D</b>		
<b>K6-344D</b>		
<b>K6-345D</b>		

STANDARD MODEL	29
BIO-BARRIER MODEL	63
Polymer's PPSU Lids	97
WINDOW MODEL	119
B <sub>3</sub> Steri MODEL	139
PLASMA MODEL	173
A1 MODEL	185
Dental, Mini, Special	205
BASKETS	223
Autolock Mechanism Emulating Indicator	243
ACCESSORIES	251

BASKET MADE FROM STAINLESS STEEL WIRE MESH

WIRE BASKET

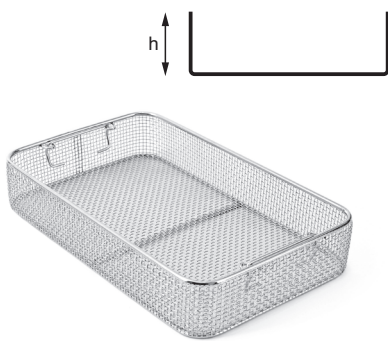

code	outside dimension
<b>K2-1500</b>	405 x 245 x h:30 mm
<b>K2-1501</b>	405 x 245 x h:50 mm
<b>K2-1502</b>	405 x 245 x h:70 mm
<b>K2-1503</b>	405 x 245 x h:100 mm



WIRE BASKET WITH LID



code	outside dimension
<b>K3-1030</b>	405 x 245 h:35 mm
<b>K3-1031</b>	405 x 245 h:55 mm
<b>K3-1032</b>	405 x 245 h:75 mm
<b>K3-1033</b>	405 x 245 h:105 mm
<b>K3-1034</b>	ONLY LID, for baskets




WIRE BASKET WITH PERFORATED METAL PLATE

WIRE BASKET



code	outside dimension
<b>K2-1501 P</b>	405 x 246 x h:50 mm
<b>K2-1502 P</b>	405 x 246 x h:70 mm
<b>K2-1503 P</b>	405 x 246 x h:100 mm

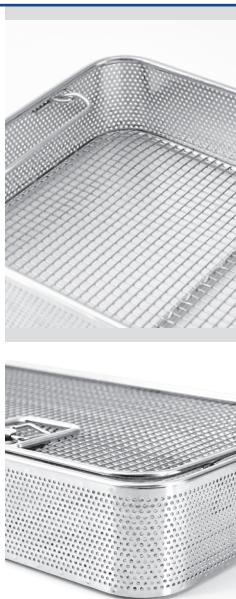

WIRE BASKET WITH LID







code	outside dimension
<b>K3-1031 P</b>	405 x 246 x h:55 mm
<b>K3-1032 P</b>	405 x 246 x h:75 mm
<b>K3-1033 P</b>	405 x 246 x h:105 mm
<b>K3-1034 P</b>	ONLY LID, for baskets






BASKET MADE FROM STAINLESS STEEL SHEET

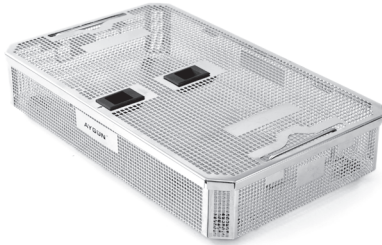BASKET (with stand)

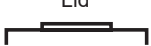


	code	outside dimension
Box and Lid 	<b>K3-728</b>	411 x 259 h:44 mm
	<b>K3-729</b>	411 x 259 h:64 mm
	<b>K3-730</b>	411 x 259 h:84 mm
	<b>K3-731</b>	411 x 259 h:114 mm
Box 	<b>K3-748</b>	405 x 253 h:40 mm
	<b>K3-749</b>	405 x 253 h:60 mm
	<b>K3-750</b>	405 x 253 h:80 mm
	<b>K3-751</b>	405 x 253 h:110 mm
Lid 	<b>K3-762</b>	ONLY LID, for baskets

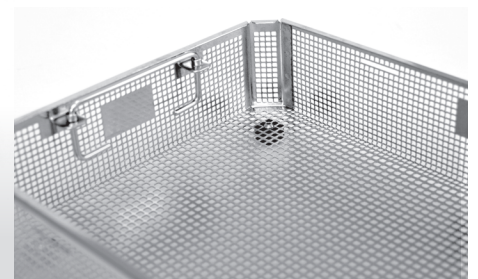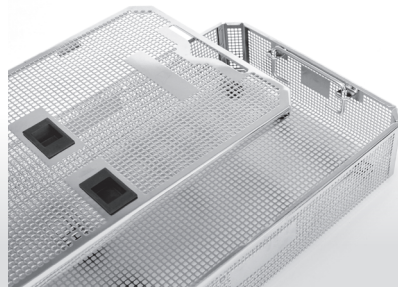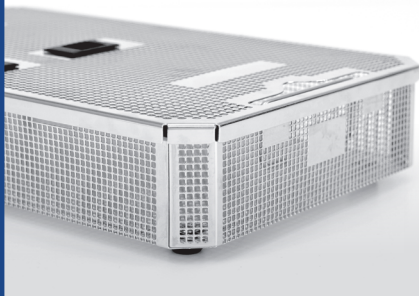

The baskets which are without lid and less than 44 mm in height do not include the handles while the ones with lid include the handles.

STAND

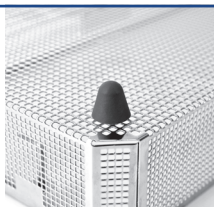

code	height
<b>K3-767</b>	20 mm

STAND



code	height
<b>K3-768</b>	10 mm

STAND



code	height
<b>K3-766</b>	5 mm

IDENTIFICATION (aluminium)

LABEL HOLDER



code	color
<b>K2-560</b>	Red
<b>K2-561</b>	Blue
<b>K2-562</b>	Green
<b>K2-563</b>	Yellow
<b>K2-564</b>	Grey
<b>K2-565</b>	Black

LABELS



code	color
<b>K2-570</b>	Red
<b>K2-571</b>	Blue
<b>K2-572</b>	Green
<b>K2-573</b>	Yellow
<b>K2-574</b>	Grey
<b>K2-575</b>	Black

Coloured identification labels  
Engraved to your specification



Grey coloured recommended for use in Sterrad.  
Color fading By effect of Hydrogen Peroxide Gas Plasma on other colors.

STANDARD MODEL 29

BIO-BARRIER MODEL 63

Polymer's PPSU Lids 97

WINDOW MODEL 119

B<sub>3</sub> Barrier MODEL 139

PLASMA MODEL 173

A1 MODEL 185


Dental, Mini, Special 205

BASKETS 223

Autolock Mechanism Emulating Indicator 243

ACCESSORIES 251

**DIVIDERS**

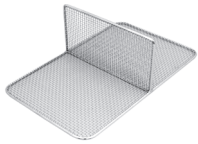


**K2-1540**  
405 x 250 x 80 mm  
**K2-1541**  
405 x 250 x 110 mm  
**K2-1542**  
405 x 250 x 130 mm

**WIRE BASE**



**K2-1522**  
405 x 255 x 30 mm


**DIVIDERS**

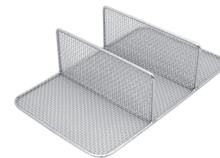

**K2-1530**  
405 x 250 x 80 mm  
**K2-1531**  
405 x 250 x 110 mm  
**K2-1532**  
405 x 250 x 130 mm

( - )

**DRAPE RETAINER**

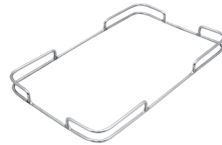

**K2-1521**  
405 x 255 x 30 mm

**DRAPE RETAINER**



**K2-1520**  
405 x 255 mm

**ALUMINIUM BASKET (with lid)**



**K3-1119**  
245 x 130 x 55 mm

**ALUMINIUM BASKET (with lid)**




**K3-1120**  
245 x 130 x 70 mm

**ALUMINIUM BASKET (with lid)**




**K3-1121**  
245 x 130 x 105 mm

**SILICONE MAT**

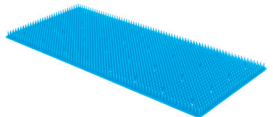

**K2-524**  
380 x 230 mm  
**K2-525**  
440 x 230 mm

**SILICONE MESH MAT**

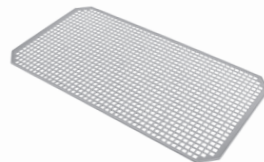

**K2-534**  
395 x 245 mm

**SILICONE HANDLE**



**K2-710**   
**K2-711**   
**K2-712**   
**K2-713**   
**K2-714**   
**K2-715**

**PAPER FILTER**




**K2-740**  
Single use 100 pcs with indicator. Aygun mark.  
• Diameter: 184 mm

**TEFLON FILTER**



**K2-742**  
Number of sterilization cycle > 2000  
• Diameter: 194 mm

**FILTER DISPENSER**

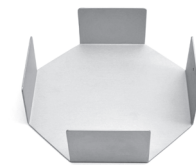


**K2-730**

**IDENTIFICATION LABELS (aluminium)**




**K2-700**   
**K2-701**   
**K2-702**   
**K2-703**   
**K2-704**   
**K2-705**

**LABELS (paper)**

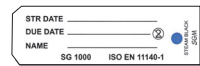

**K2-673**  
100 pcs. With indicator.  
For small size and special containers

**LABELS (paper)**

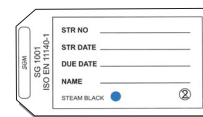

**K2-672**  
100 pcs. With indicator.  
For full and Half size containers


( - )

**AUTO LOCK MECHANISM**



**K1-900**  
(2 Pcs)

**EMULATING INDICATOR CARD (class 6)**

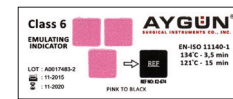

**K2-674**  
Paper 250 Pcs.



**SECURITY SEAL (blue)**  
**K2-663**  
100 pcs. with indicator. made of plastic



**SECURITY SEAL (blue)**  
**K2-664**  
100 pcs. made of plastic



**SECURITY SEAL (blue)**  
**K2-661**  
100 pcs. with number. made of plastic



**SECURITY SEAL (red)**  
**K2-662**  
100 pcs. with indicator. made of plastic



**SECURITY SEAL (red)**  
**K2-660**  
100 pcs. made of plastic



**SECURITY SEAL (red)**  
**K2-665**  
100 pcs. with number. made of plastic



**DISPOSAL STAMP**  
**K2-717**  
40 x 35 x 9 mm  
Requested number or letter must be informed during the order.



**LABEL HOLDER**  
**K2-716**  
50 x 35 x 13 mm  
Holder for identification labels

( - )